



# Conservative Treatment

Excavator	38
Filling Instrument	40
Composite Syringe	44
Amalgam Product	44
Matrix Retainer / Separator	45
Contact Protector	45
Applicator	46
Rubber Dam Product	47
Endodontic Product	50
Fiber Cutter	58





Excavator



-Round Tip-




Str. #1 / Cvd. #1

Str. #2 / Cvd. #2



Str. #3 / Cvd. #3

Str. #4L / Cvd. #4L

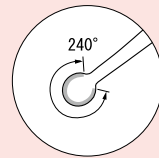\*Material / Si (Handle)




\*The edge along with the circle at a tip remove softened dentin and make cavity preparation efficiently.

Description		#1	#2	#3	#4L
Str.	 Cavity preparation for wall and floor.	01601	01602	01603	01604
Cvd.	 Used for undercut area, narrow rim of an opening and at the cavity angle.	01701	01702	01703	01704

-CDR-

\*Excavator which pursued minimal intervention.

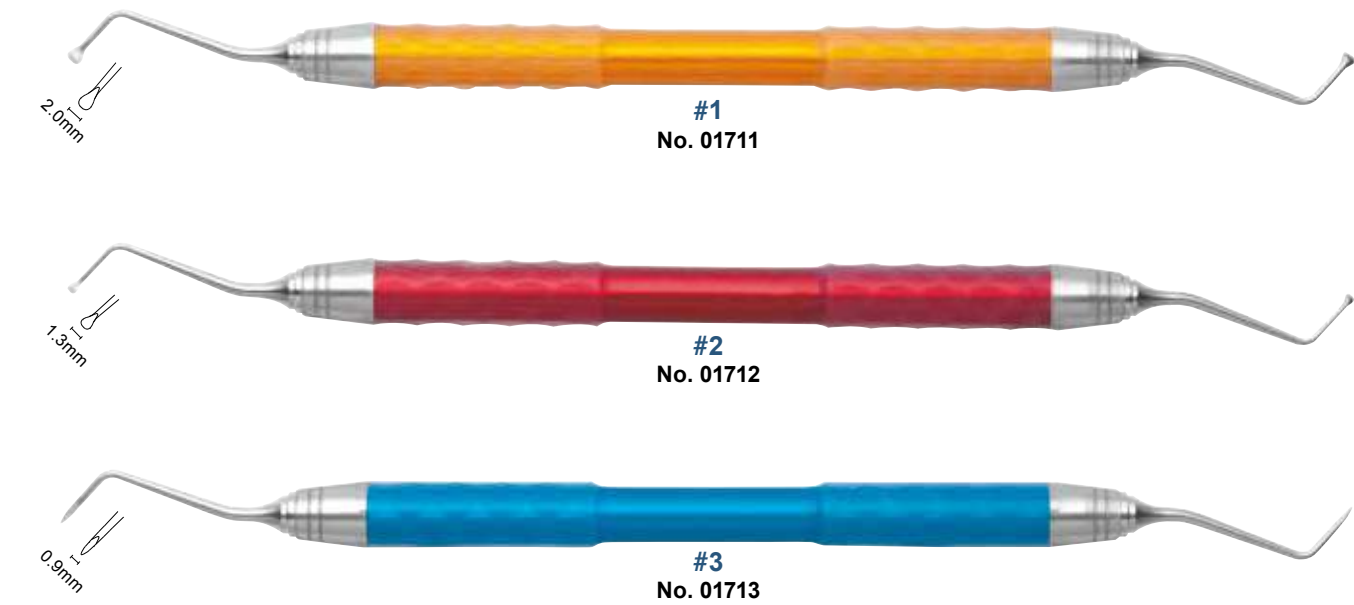


#1  
No. 01711

#2  
No. 01712

#3  
No. 01713

\*Material / Si (Handle)

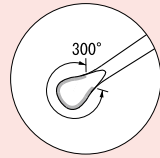

#1 / #2  
\*300° wide edge allows efficiently removal of cavity dentine.  
\*Slope working part design allows easy clean of particles.








## Composite Filling Instrument







### -TMDU 2 (Tokyo Medical Dental University)-

#### \*Diamond Like Carbon Coating

\*Remarkably stick-resistant for easy workability. Improved strength, and superior wear resistance.

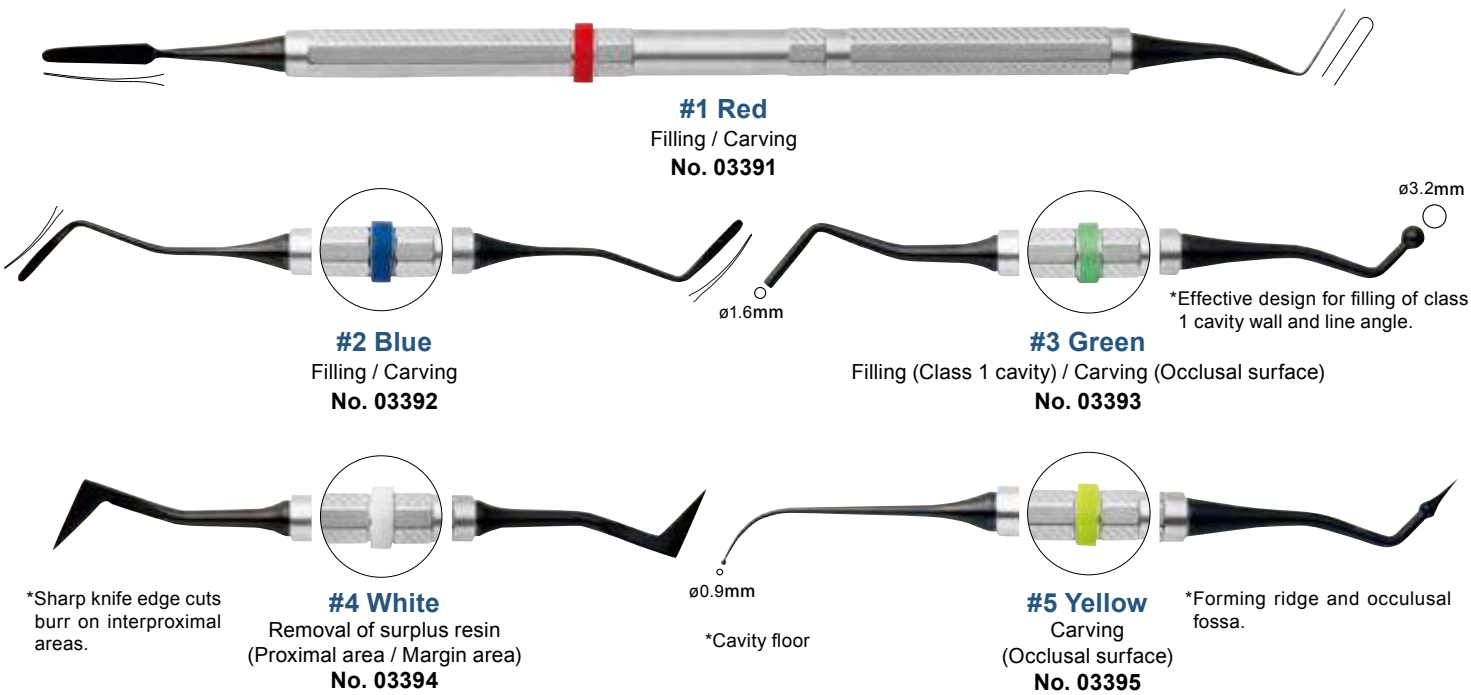





\*Thickness / 0.2mm (#1)  
0.1mm (#2)  
\*Material / Si (Color Ring)

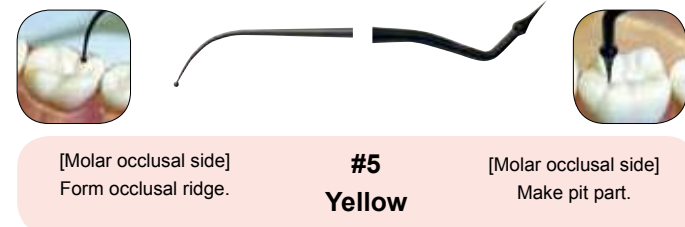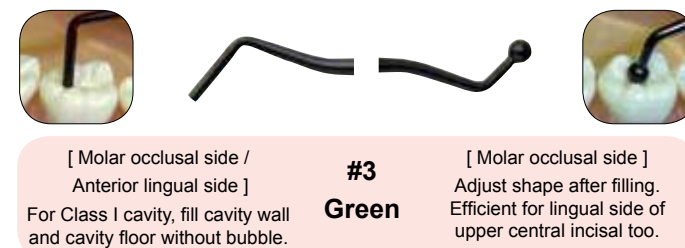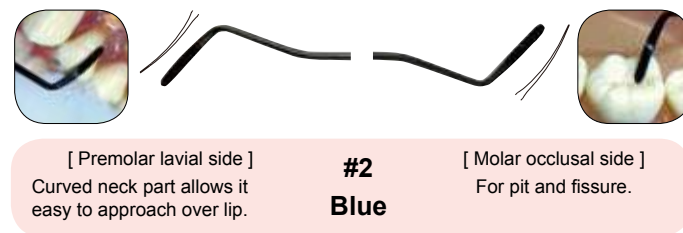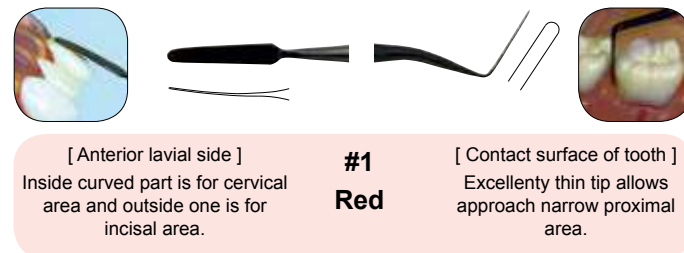

### -TMDU 3 (Tokyo Medical Dental University)-


#### \*Diamond Like Carbon Coating

\*Remarkably stick-resistant for easy workability. Improved strength, and superior wear resistance.

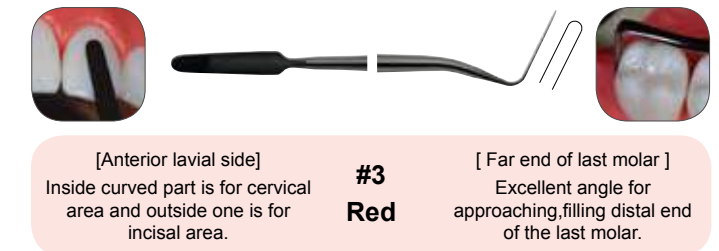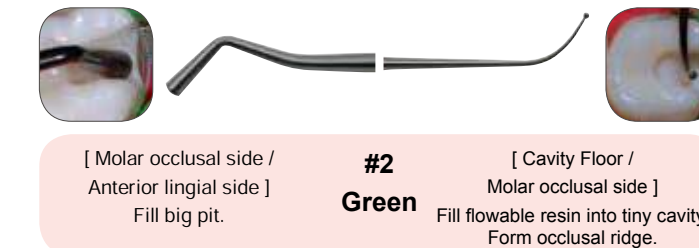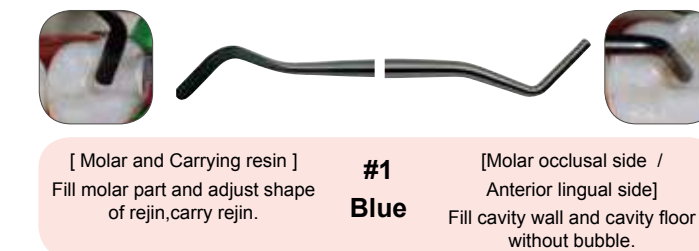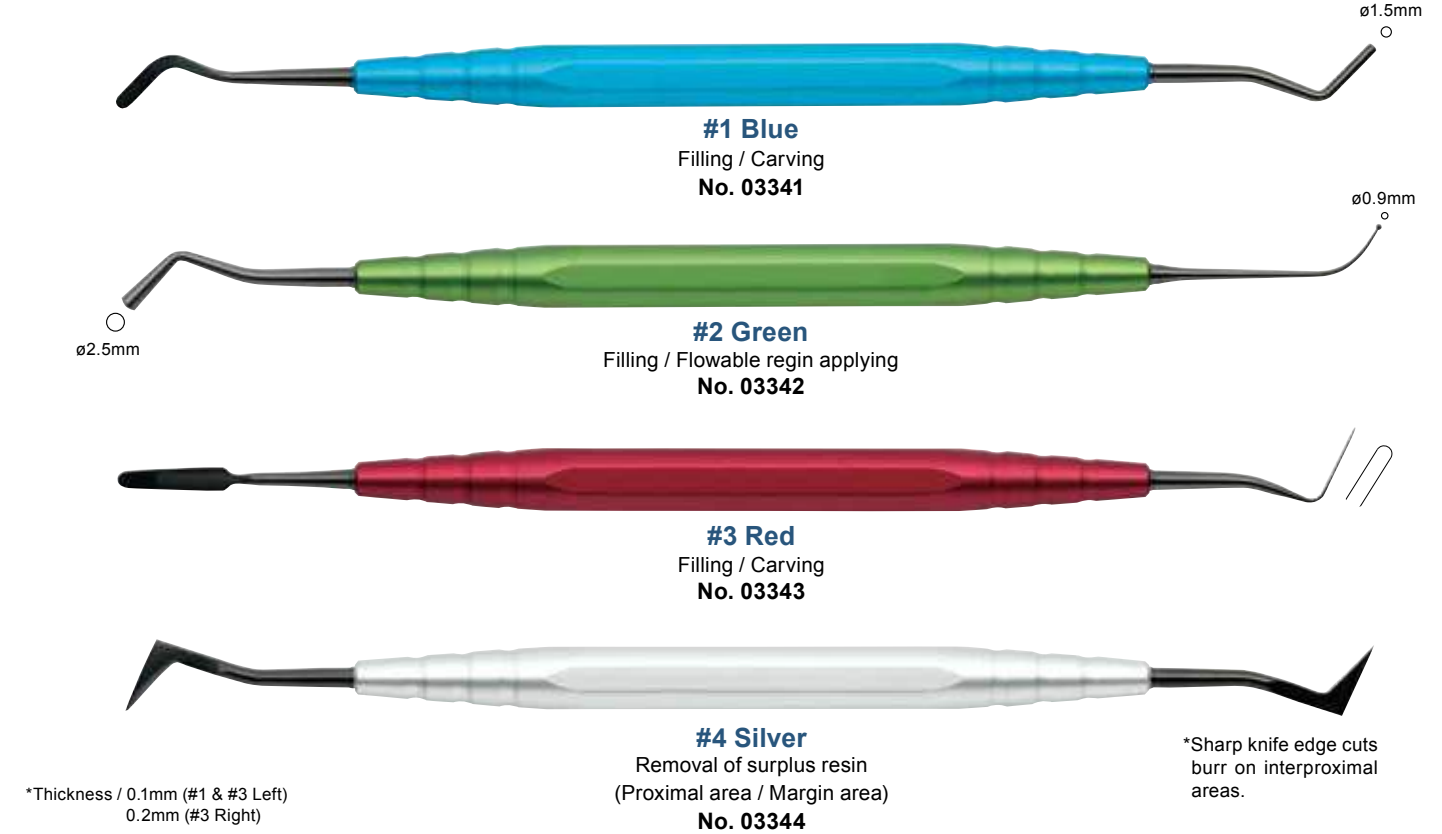

Contents

L : Main body L / Nozzle & Plug L (100 ea.)  
S : Main body S / Nozzle & Plug S (100 ea.)

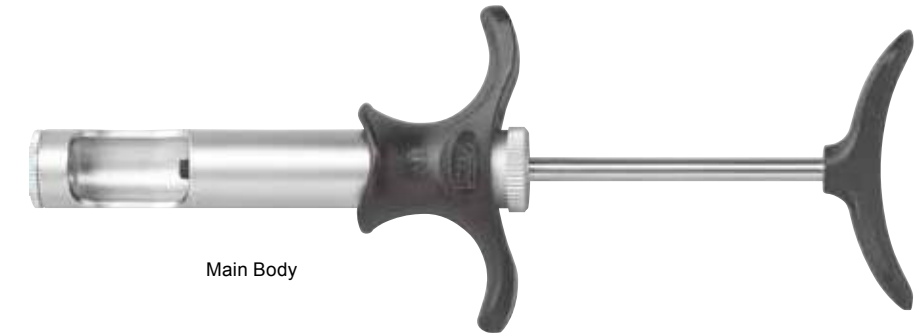

Main Body


Nozzle & Plug L  
(Plug / Blue)



L

No. 13251


Nozzle & Plug S  
(Plug / Green)


S

No. 13252

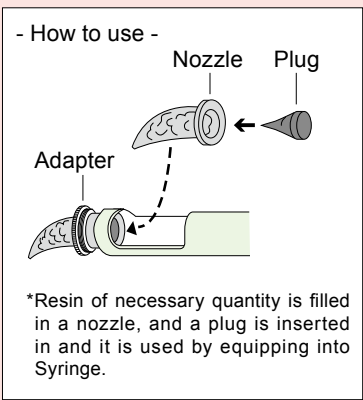

Individual Parts

No.	Description
13255	Main Body L
13256	Main Body S
13188	Nozzle & Plug L (100 ea.)
13190	Nozzle & Plug S (100 ea.)
13193	Nozzle & Plug L (500 ea.)
13194	Nozzle & Plug S (500 ea.)

\*Material / PPS (Black part)

## Amalgam Burnisher



S

No. 05002

## Amalgam Carver



Round Handle

No. 07100


## Matrix Retainer Set



Contents

Main body / Band (10 ea.)

No. 13170

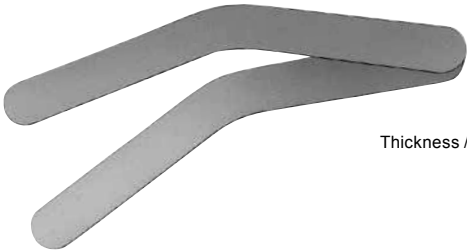

Thickness / 30μm

Individual Parts

No.	Description
13171	Main Body
13172	Band (10 ea.)

- How to use -


1) The both ends of a matrix retainer band are united and made a ring.

Crown part




Dental neck part

2) The matrix retainer band which became a ring is inserted in a matrix retainer and a fixed screw is tightened.

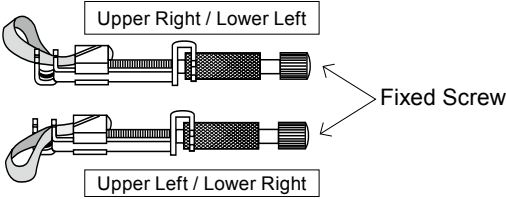

3) Apply to the target tooth and tightened adjustable screw.

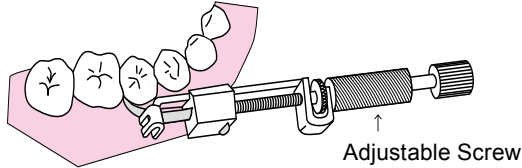

## Separator



Ivory

Anterior

No. 13161



Elliot

Posterior

No. 13162



## Contact Protector Set

Contents

Anterior (10 ea.) / Posterior (10 ea.)

\*Stainless Steel


\*Protect adjacent tooth during the cavity cutting.

No. 13175

Anterior




Posterior

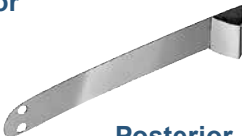

\*Thickness / 50μm

\*Not sold separately.

- How to use -

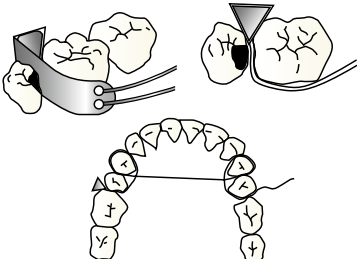

\*It is used after letting froth pass in a hole.




## Lining Applicator



**-Round-**



**-Octagonal-**

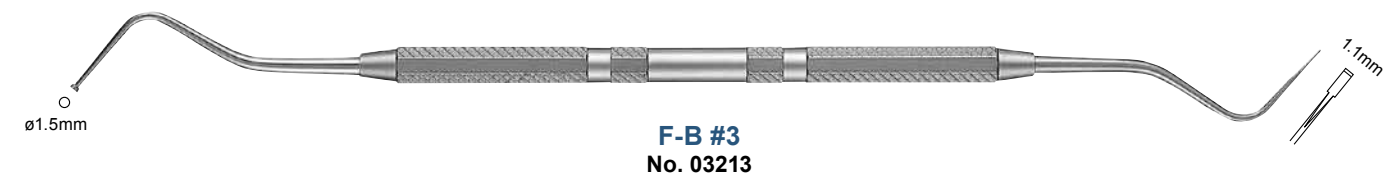

## Ball-Pointed Applicator



No. 25015





## Rubber Dam Clamp




**\*Stainless Steel**

\*#12A / 13A → Serrated jaws and offset jaw height provides stability to partially erupted teeth.

\*#26 / #26A → Used in pedodontic. In that case, #26 is for deciduous molars and #26A is for deciduous anterior teeth.

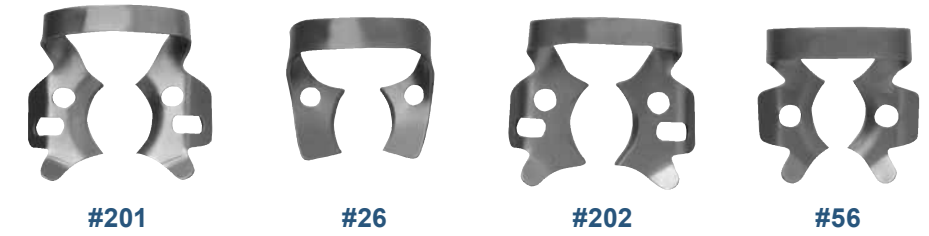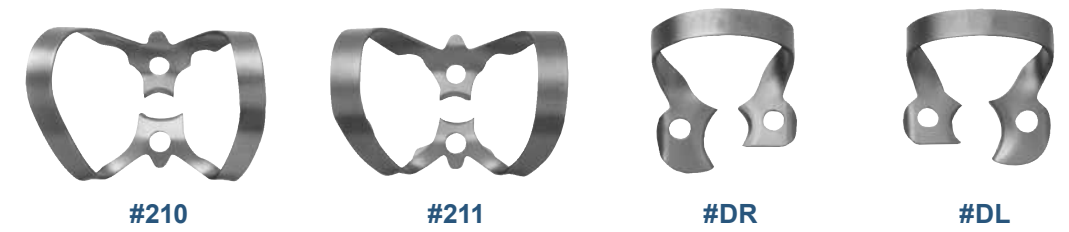

-Clamp shape-

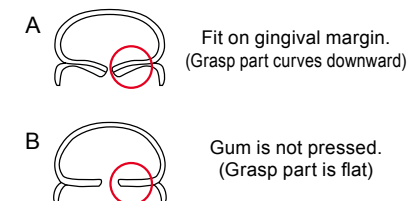

No.	Description		Shape
08526	#26	Upper molars	A
08501	#201	Upper molars	A
08556	#56	Lower molars	A
08502	#202	Lower molars	B
08506	#206	Premolars	A
08507	#207	Premolars	B
08527	#26A	Premolars	B

**Board  
No. 08258**



Size / W132 x D87mm

- How to use -








## Reamer Guard

- \*Can be arranged "Reamer Guard" by reamer or file type and size.
- \*Cover will not allow the reamers and files to fall out.
- \*An Engine reamer can be stored with "Bur Guard".

No.	Description
08431	Yellow
08432	Orange
08433	Pink
08434	Blue
08435	Green



### -Holding length-

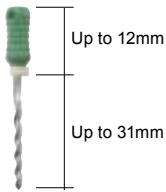

- \*Storage capacity / 8 reamers and files
- \*Size / W74 x D10 x H52 mm
- \*Material / POM
- \*Sterilizable with liquid sterilizing solutions.
- \*Product does not include reamers and files.

### Sold Separately


#### Case No. 08439



- \*Storage capacity / 18 ea.
- \*Size / W167 x D126 x H62mm
- \*Sterilizable with liquid sterilizing solutions.
- \*Product does not include "Reamer Guard".



\*The cover becomes a stand.

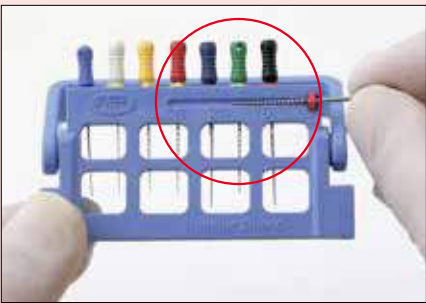

\*Has a built in gauge to measure stop positions during root canal treatments.




\*Able to connect each "Reamer Guards".

## Reamer Stand

### -Carry Set-

**Contents** Carry (3 ea.) / Case

#### No. 08411

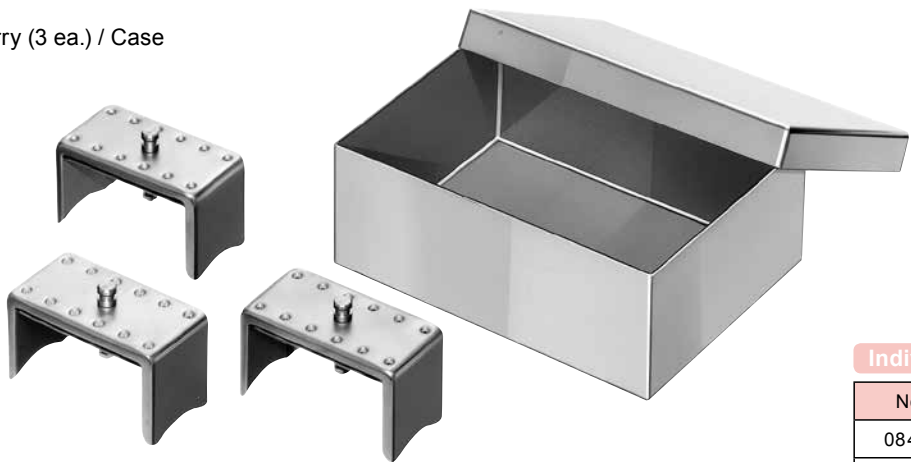

- \*Storage capacity / 12 reamers and files (Carry)  
/ 6 ea. (Case)
- \*Size / W60 x D30 x H33mm (Carry)  
/ W120 x D90 x H40mm (Case)


### Individual Parts

No.	Description
08412	Carry (3 ea.)
24016	Case

### -Safety-

- \*It is ideal for managing reamers and files for every patients.
- \*Can be autoclave or ultrasonic cleaning with reamers and files.
- \*Lid does not allow reamers and files to fall out.

#### No. 08415

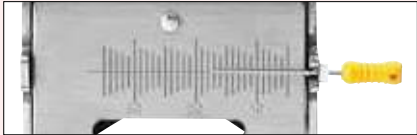

\*Has a build in ruler.

- \*Storage capacity / 11 reamers and files
- \*Holding length / From 21 to 31mm
- \*Size / W48 x D43 x H47mm
- \*Material / PP (Green part)
- \*Product does not include reamers and files.



## Clean Pot

**Contents** Main Body / Paper (100 sheets)

#### No. 08423

- \*During the treatment, can be cleaned reamers and files by paper.
- \*Set 2 sheets to stabilize reamers for cleaning.



### Individual Parts

#### Paper No. 08422 (300 sheets)




- \*Size / ø45 x H36mm
- \*Product does not include reamers and files.






Contents Main Body / Joint (ø10mm & ø11mm Tapered) / Nozzle (5 ea.) / Cleaning Wire

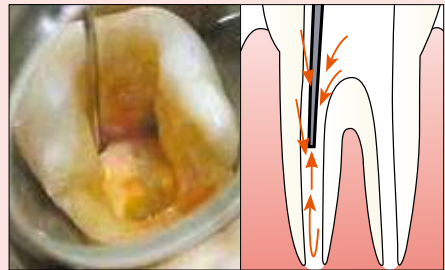

\*The inside of a root canal is dried by an air current with suction in a root canal.  
\*Because of drying by an air current, it also becomes saving of a paper point.

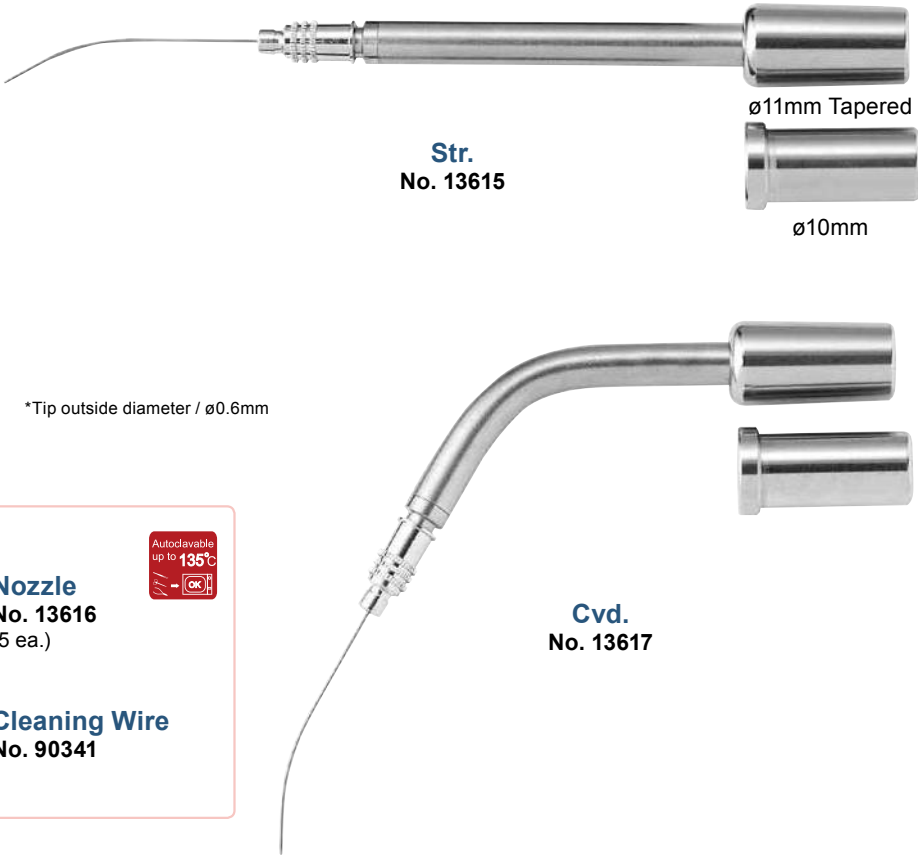

Individual Parts



Sold Separately

Cleaning Brush S (ø4mm)  
No. 13651  
(12 ea.)



Adapter  
No. 13614



\*In case the inside diameter of evacuation tube is 16mm, this adaptor is needed for the use of Evacuation Tip Universal.


\*Material / PTFE

## Irrigating Syringe



\*Material / PFA (Cap / Main body)  
PP (Trigger)

## Explorer

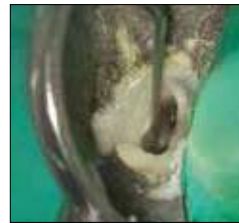

## Broach Holder

No. 13150



## File Holder Type-2

\*Excellent for use in hard-to-reach posterior teeth.  
\*Applicable for root canal treatment of both mesial and distal.

No. 08064



## Finger Ruler

No. 21155




## Scale



## Endodontic Gauge

No. 21154

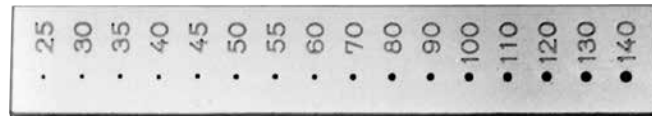



Front side



Back side





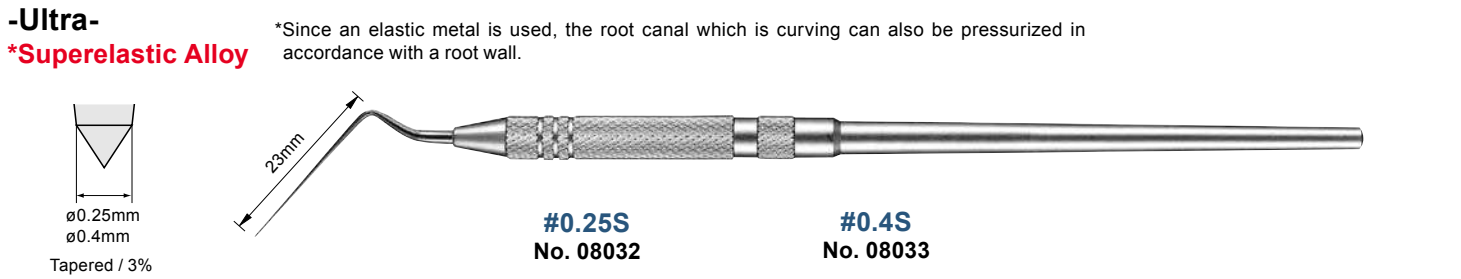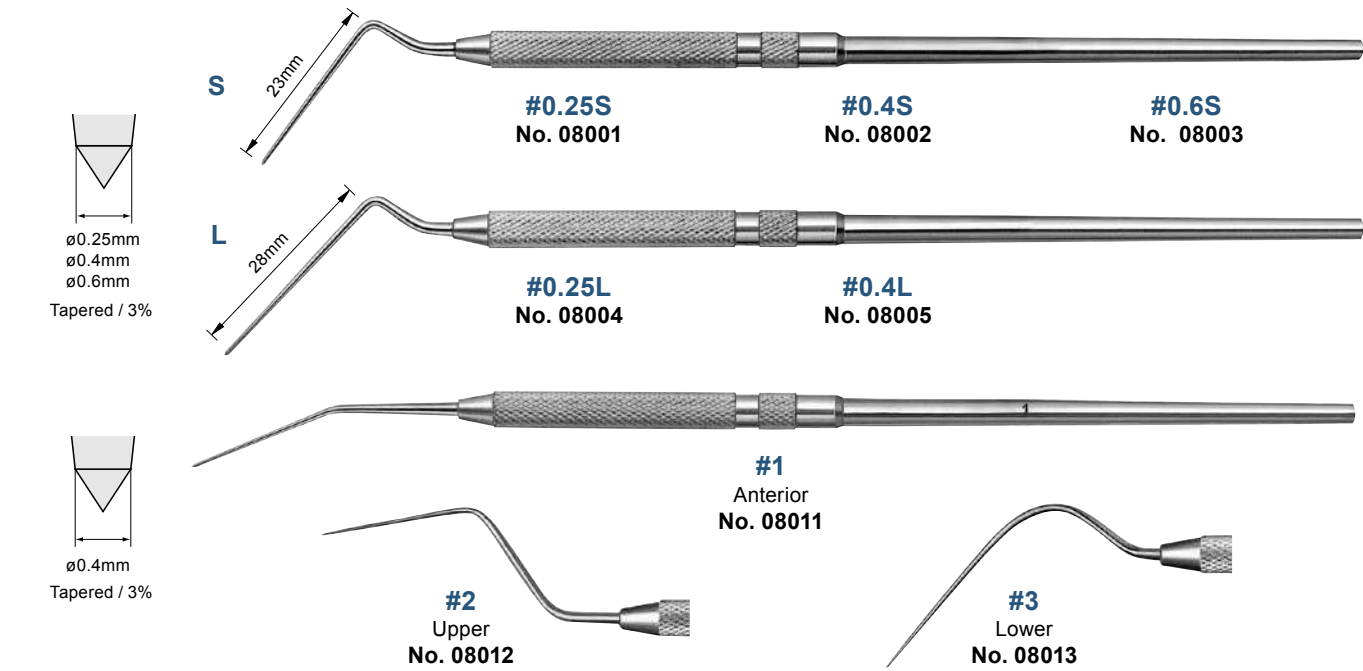


## Spreader Tip (Replaceable Type)

\*Remove all joint parts and replaceable tip before sterilization and disinfection.  
\*For more replaceable type products see P. 201.



### -Handle-

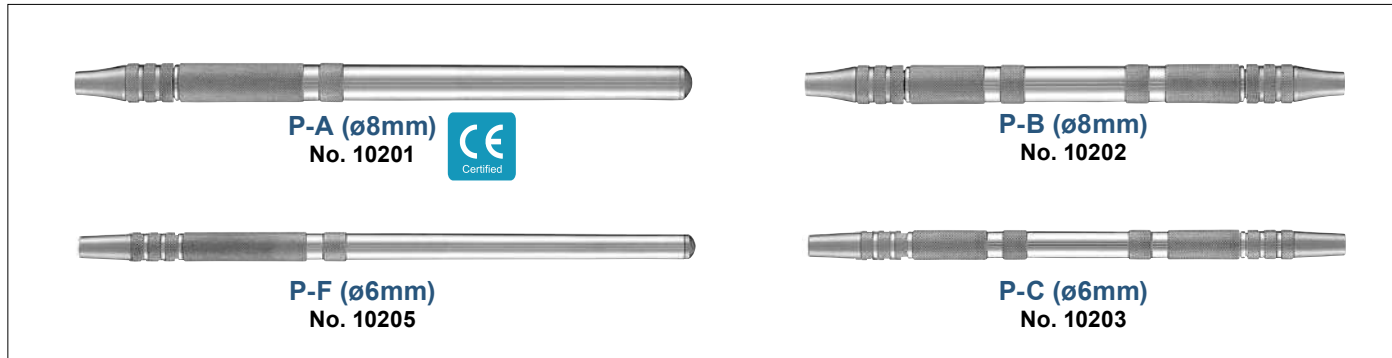


### -Tip (3 ea.)-




## Root Canal Plugger

### -Round-

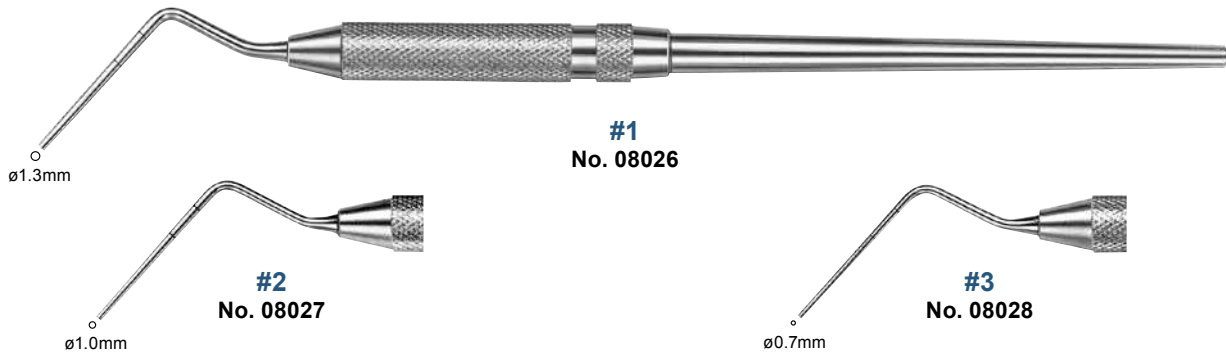


### -Octagonal-

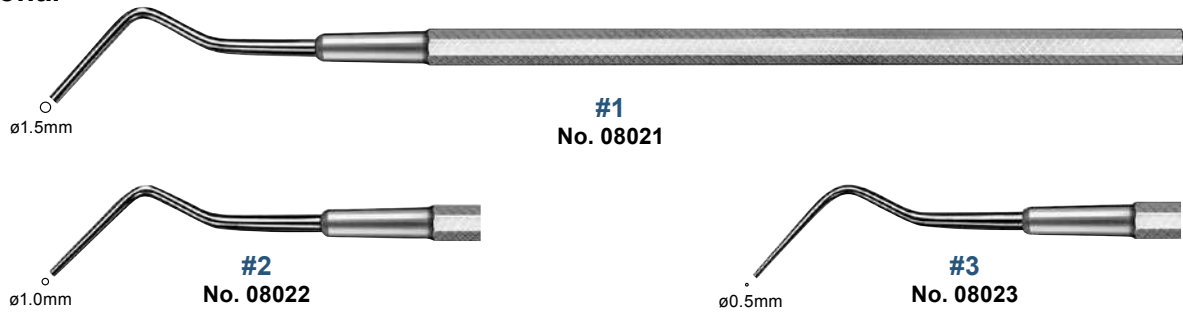


### -with Heat Cutter-

\*Combination of heat cutter and plugger makes cutting and compacting of extra gutta percha easily.  
\*Heat cutter is designed with a metal heat retention bulb for longer working time.

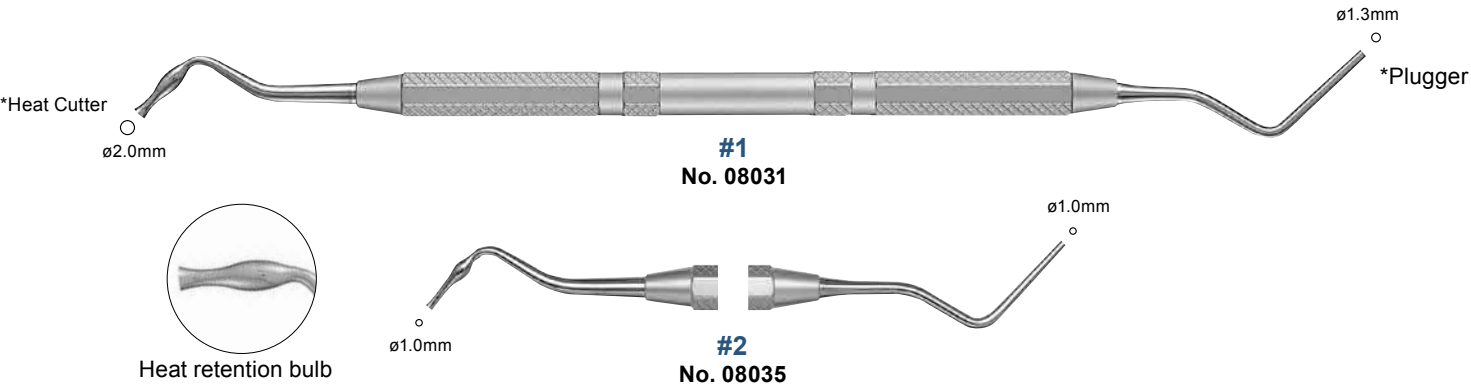

## Heat Cutter

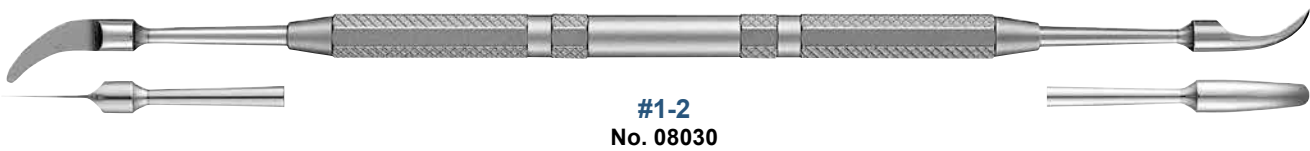





## Fiber Cutter Type2

**Contents** Main body / Holder / Sleeve Guide (3 ea.)

- \*Easy to cut the fiber post length for your need.
- \*There is no fluff of a cutting plane.
- \*Particulates are not generated.





**No. 08071**



- \*Cut fiber can be entered in to the holder.
- \*Easy to cut the fiber thanks to the gauge on the holder.  
(Measurement from 5mm to 15mm)

\*Size / 130mm  
\*Material / TPX (Holder)

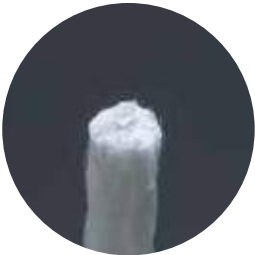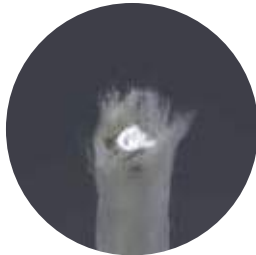

Fiber Cutter Type2	Nipper	Steal Round Bur	Cutting Disc	Cutting Disc with Diamond
				

### Problems from using Steal Round Bur, Cutting Disc, Diamond Disc

- 1) Fiber become fluffy.
- 2) Melt the surface because of the friction heat between fiber post and bur.
- 3) Remain the color of bur or disc on the section.